

SPACE SAVING CHART

SSC-350/350CG

Eye & Health Care




SSC-350 / 350CG Specifications

Refracting Distance	5 m
Installation Distance	1.1 m (from the chart window of SSC-350 to patient's eye)
Charts Changeover Devices	Remote control (Option) or control box of refractor
Chart Switching / Masking Time	0.5 seconds (average)
Masking Device	Horizontal line, Vertical line, Single letter
Size of Chart Window	338 (W) x 258 (H) mm
Lamp	12 V, 30 W (Halogen lamp)
Dimensions / Weight	Main body: 396 (W) x 270 (D) x 560 (H) mm / 18 kg, 20 kg (CG) 15.6 (W) x 10.6 (D) x 22 (H) " / 40 lbs., 44 lbs. (CG) Remote controller: 66 (W) x 217 (D) x 19 (H) mm / 198 g 2.6 (W) x 8.5 (D) x 0.7(H) " / 0.4 lbs.
Power Source	AC / 100 / 120 / 220 / 230 V, 50 / 60 Hz, 110 VA, 130 VA (CG)

1. The SSC-350 requires only 1.1 m to test at a 5 m distance. Limited space is not a problem anymore.
- *2. The ratio of the contrast is adjustable in 4 levels : 100%, 25%, 12.5% and 6%.
- *3. The contrast measurement can be carried out under the high intensity of the white LED glare lamp.
4. Chart for high visual acuity : 2.5.
5. The SSC-350 incorporates the Log MAR (Logarithm of the Minimum Angle of Resolution) chart which is highly effective for follow-up of patients after refractive surgery.

*SSC-350CG only

26 Built-in Charts

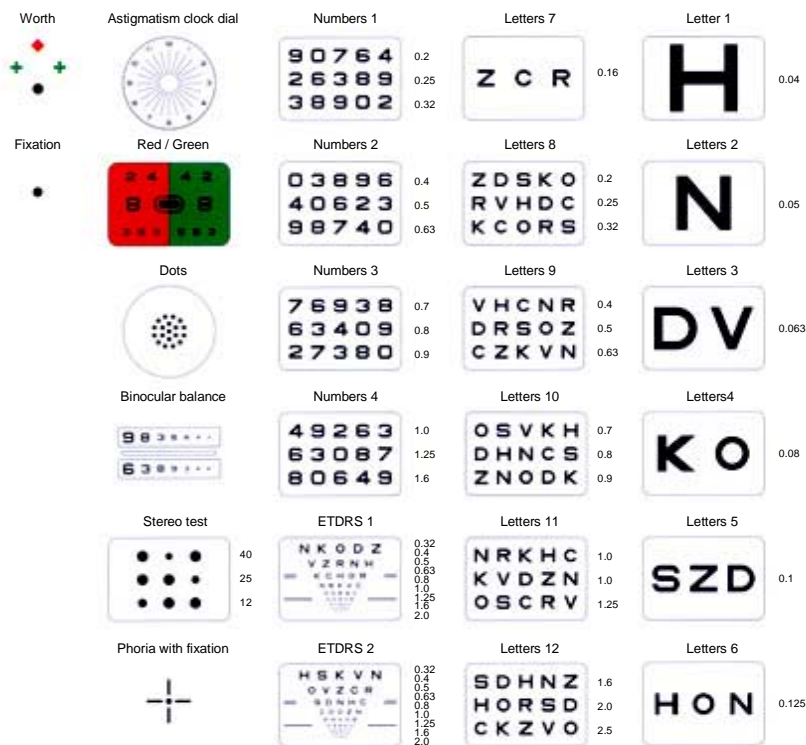

Specifications and design are subject to change without notice for improvement.



SSC-350



NIDEK CO., LTD.

HEAD OFFICE

34-14, Maehama, Hiroishi
Gamagori, Aichi 443-0038, Japan
Telephone :81-533-67-6611
Facsimile :81-533-67-6610
URL :http://www.nidek.co.jp

TOKYO OFFICE

(International Div.)
6F Takahashi Bldg.,
3-2 Kanda-Jinboucho
Chiyoda, Tokyo 101-0051, Japan
Telephone :81-3-3288-0571
Facsimile :81-3-3288-0570
URL :http://www.nidek.com

NIDEK INC.

47651 Washington Drive
Fremont, CA 94539, U.S.A.
Telephone :1-510-226-5700
:1-800-223-9044 (US only)
Facsimile :1-510-226-5750
URL :http://www.usa.nidek.com

NIDEK TECHNOLOGIES AMERICA INC.

5500 West Friendly Ave.
Suite 101 Greensboro, NC 27410, U.S.A.
Telephone :1-336-851-0225
:1-888-382-5064 (US only)
Facsimile :1-336-851-0917
URL :http://www.nidektech.com

NIDEK SOCIETE ANONYME

Europarc
13, rue Auguste Perret
94042 Creteil, France
Telephone :33-1-49 80 97 97
Facsimile :33-1-49 80 32 08
URL :http://www.nidek.fr

NIDEK TECHNOLOGIES SRL

Via dell'Artigianato, 6/A
35020 Albinasego (Padova), Italy
Telephone :39 049 8629200 / 8626399
Facsimile :39 049 8626824
URL :http://www.nidektechnologies.it